


Digital Microscope and
Measurement System

INSPEX 3

Smart Inspection

ASH

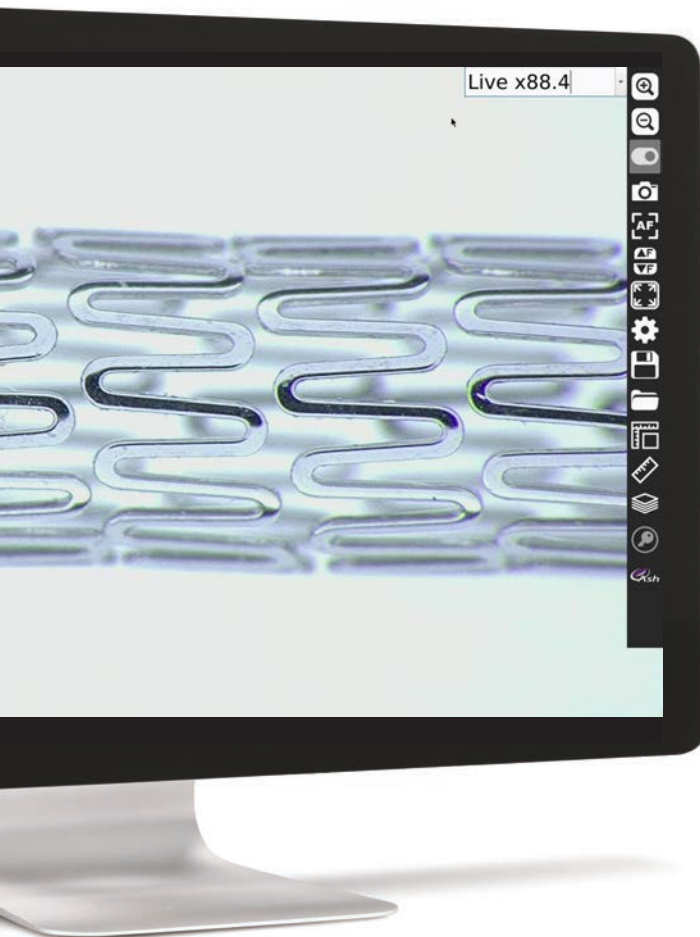

INSPEX 3

Smart Inspection



Simply Smarter

INSPEX 3

Digital Microscope System


Inspex 3 is a powerful, flexible and intuitive HD digital microscope designed around ASH's new 30X HD Camera module.

This zoom camera module was designed by ASH to meet industrial inspection needs.

ASH has developed unique patented high-speed image processing algorithms deliver exceptional image quality, super-fast focus speed, and smooth digital zoom that makes inspection simpler, ergonomic and efficient.

Inspex 3 is ideal for industries such as medical device, pharma, precision engineering and electronics.

The benefits of the Inspex 3 system are clear to see.

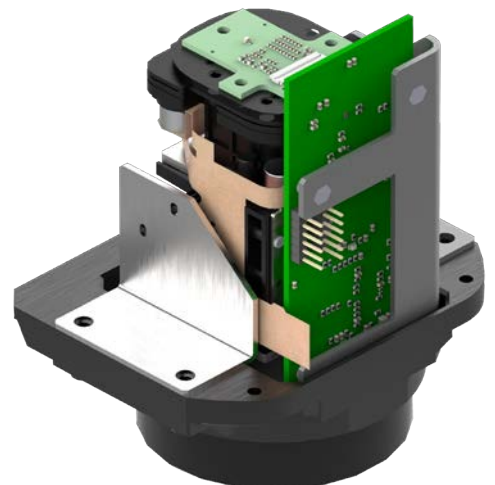

INSPEX 3

Smart Inspection

New Features

Inspex 3 is the latest evolution from ASH. Incorporating the AshCam™ module with a wide range of additional and improved capabilities.

Superb Image Quality


Experience unrivalled Full HD video imaging never seen before in an Ash system.

The Inspex 3 offers an enhanced, vibrant, crystal clear image for even the most demanding inspection applications.



Super Fast Auto-Focus

Place any part under the camera and it will immediately adjust the focus position allowing you to quickly and seamlessly inspect your part without having to adjust manual focus or change the height of the part.

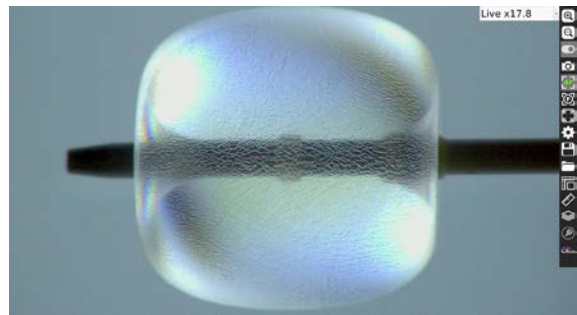

Manual Rocker Focus

Save time by quickly adjusting the focus using the rocker icon to inspect the area of interest on a sample part when in manual focus mode.



SpotFocus™

Quickly focus on the area of interest by simply using the mouse pointer.

Improve accuracy and reduce human error with the capabilities of Spot Focus.*



* Global patent applied for by ASH.


AshTruColour™ - True Colour Reproduction

View true colour reproduction of your sample with the Inspex 3. Replicate real and accurate colours as seen with the naked eye for true colour representation.




Advanced Camera Settings

The new Advanced Camera Settings gives the user more power to enhance the image for a wide range of inspection and measurement capabilities. Tailor sharpness, contrast, saturation and camera shutter speed to suit your particular application.

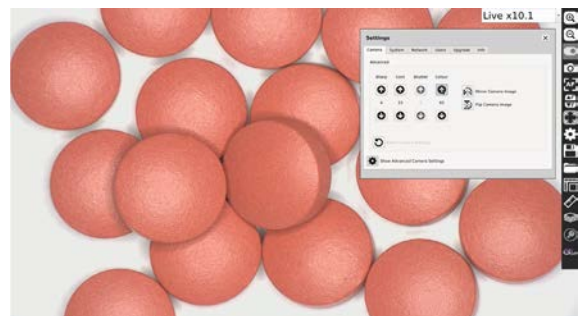


RTLDC™

Real-Time Lens Distortion Correction*. Lens distortion is inherent in all microscopes. Image distortion at the outer edges of large samples is automatically corrected by the Inspex 3.

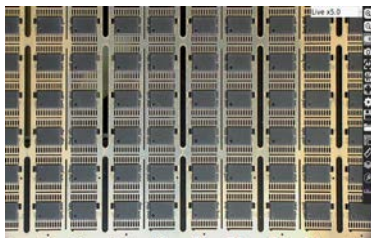


Inspex 3 RTLDC

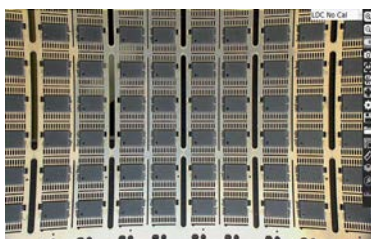

Without RTLDC

No Video Latency

The video from the Inspex 3 has zero delay. It's 3x times faster than our previous systems. View parts in real time with no video lag, allowing you to comfortably inspect, rework, modify and assemble any part. There is no delay between movement under the camera and what you see on the screen.

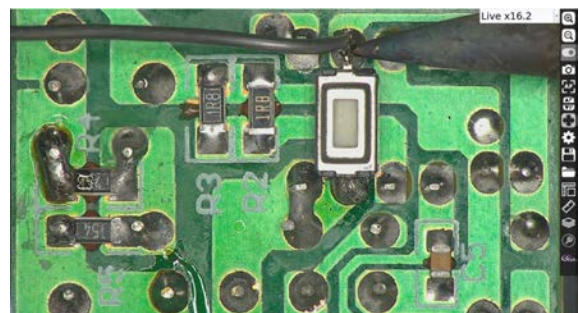

* Global Patent applied for by ASH.

INSPEX 3

The Inspex 3 maintains the following improved great features introduced with the hugely successful Inspex II.



AshCal™

Inspex 3 is factory calibrated before shipping. No time is wasted performing recalibrations between changing magnifications.



2D Measurement, Annotation & Graticules

Point to point measurement and annotation of samples and creation of graticules.



User Privileges

User privilege settings enables operational control and traceability. Assign multiple users with access to different settings and features, improving security and streamlining the inspection process.



Interchangeable Lighting

Multiple types of interchangeable lighting available, including Ring Light, DomeLight, Polarised Light and UV Light to address any inspection application.



Image Stacking

Inspex 3 automatically captures several images at different focal depths to create a fully focused, sharp, clear image for easy inspection. Save time by removing the need to adjust the camera height or manual focus.



Save to Network

Networking enables direct saving to the server or cloud for increased workflow efficiency.



On Screen Preset Buttons

On Screen Preset Buttons allow quick access to pre-configured part-specific system settings.



Image Stamping

Image Stamping with time, date, user and magnification level. Easy documentation & traceability for accurate quality control records using image capture to USB key.

INSPEX 3

is the latest evolution

AshCal™

Automatic Calibration
Tracking

Improved

Auto-Focus Speed

New! Spot Focus

Rocker Focus

No Video
Latency

Improved
Depth of Field

Improved
Sensitivity

Smooth
Magnification
Zoom

RTLDC™



Simply Smarter



PLUGIN
AND **GO!**



INSPEX 3

Digital Microscope System

Smart Inspection

Applications

- In Process Visual Inspection, Rework and Assembly of Component Parts
- Connected for Traceability and Documentation. Store images and data on the Cloud.
- Large Scale Inspection in a Production Environment
- Final QC Release Inspection

Industries

- Medical Device Manufacture
- Pharma
- Precision Engineering
- Electronics

Steven Martin,
Dakin Flathers

Dakin Flathers

"The ease of use has meant that we have had no issues with operator engagement and everyone who we asked insists it has taken away any degree of uncertainty that was always evident in analogue checks.

From a personal point of view, right from the initial demonstration through to the after-sales service, I cannot fault the staff at ASH, who have gone over-and-above to assist however they have been able. Possibly the greatest testament to the unit is that other sections within the factory have seen the potential it brings and have begun making their own case to me for the purchase of additional units. I have to say, it is only a matter of time!"

Pádraig McVeigh, Creganna
- TE Connectivity

Conductix Wampfler

"ASH supplied a digital microscope system to help assessment of tube finishes without major impact to the overall cycle of production. ASH have provided excellent support from problem classification, through vision system development, right to fixturing and set-up support.

The speed of the response and ability to meet our tight timeline was very impressive.

The availability and assistance from team members of all levels, from Sales to Engineering, was invaluable to our team!"

Roland Rücker, Quality
Engineer Optoelectronics,
Heraeus Noblelight

"Our goal was to acquire a microscope with the smallest possible surface differences between 5 and 500 microns without fatigue. Likewise, we wanted to create sharp images with the abnormalities and these surveying. The implementation of these requirements has now been confirmed in practice by our employees. What also stands out is the intuitive user interface our employees were able to use within a very short time. The support from their side for commissioning and technical questions is also to be positively emphasized. We can highly recommend this microscope because of our positive experiences and you as a partner."





INSPEX 3

FI 806-001

System Components Included



LENSES

+5 Lens



LIGHTS

LED Ring Light



OTHERS

Wireless Keyboard and
Mouse

USB Memory Stick

HDMI Cable

Optional System Components



LENSES

+5 Lens

AI 280-150

+10 Plan 1x Lens

AI 100-055

+25 Lens

AI 100-053

Ash 360 Rotating Viewer

AI 801-422

Polarising Lens
(Analyser)

AI 100-041

Sub-Stage Polarising
Film & Analyser

AI 801-835

Sub-Stage Polarising
Film

AI 801-836

Polarised Ringlight &
Analyser (58mm)

AI 801-423



LIGHTS

Diffused LED Dome Light

AI 100-045

UV Ring Light 367nm

AI 801-421

Polarised Ringlight &
Analyser (58mm)

AI 801-423



CONTROLLERS

KPII External Keypad

FI 806-002

KIII External Keypad

FI 806-003



OTHERS

Ash PC Capture

AI 100-052

24" Monitor

AI 801-416

Hard Carry Case

AI 801-563



STANDS & STAGES

XY Stage for
Uplight

AI 100-011

XY Stage

AI 100-010

Large XY Stage

AI 100-057

Illuminated Track
Stand

AI 100-036

Track Stand

AI 100-037

Articulated Arm
Stand


AI 100-039

Dual Arm Boom
Stand

AI 100-038

INSPEX 3

Digital Microscope System

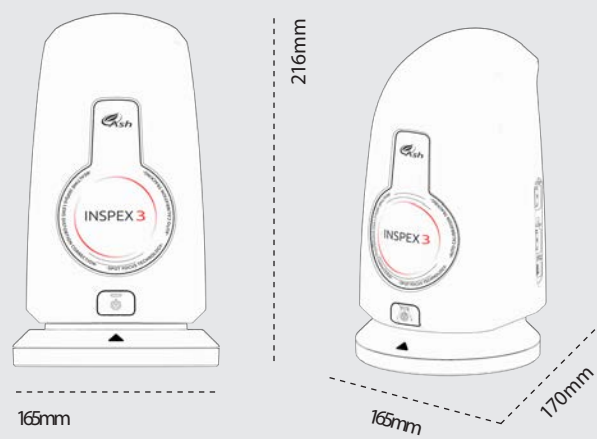

Magnification

Lens Type		+5	+10 Plan 1x	+25	+50
Optical	Magnification Range (X)	2.1 - 65.6	4.1 - 130	51.7 - 323.4	97.7 - 628
	X-axis FOV (mm)	240 - 8.34	80 - 4.1	3.5 - 1.75	2.1 - 1
	Y-axis FOV (mm)	135 - 4.69	71 - 2.2	1.75 - 0.9	2 - 0.6
Digital	Magnification Range (X)	66.1d - 131.6d	130.5d - 259.9d	325d - 646.7d	631.2d - 1256.1d
	X-axis FOV (mm)	8.34 - 4.13	4.15 - 2.08	1.65 - 0.85	1.2 - 0.44
	Y-axis FOV (mm)	4.69 - 2.32	2.32 - 1.16	0.9 - 0.5	0.7 - 0.26
	Working Distance (mm)	195	79	49	36
	Depth of Field (mm)	80 - 0.5	42 - 0.2	0.3 - 0.015	15um - 4um
	Video Latency (milli seconds)	20 (50Hz)/ 17 (60Hz)	20 (50Hz)/ 17 (60Hz)	20 (50Hz)/ 17 (60Hz)	20 (50Hz)/ 17 (60Hz)

Technical Specifications

Zoom Range (with supplied +5 Lens)	INSPEX 3
Camera Resolution	2.1 - 131.6x
Monitor Connections	1920 x 1080 pixels
Monitor Requirements	HDMI / DVI
Input / Output	HD Ready / Full HD (Recommended)
	HDMI Output
	USB 2.0 (x4 Ports)
	Mini USB Port
	General Purpose IO (x3 Ports)
	DC Power Jack 24V
Internal Storage	16GB
Image Capture	Internal Storage
	Removable USB Image Storage
	USB on the Go (PC Connectivity)
Power	24W
Dimensions	216mm x 165mm x 170mm
Weight	1.5kg
Operating Temperature	Storage -10°C to +60°C
	Operating +5°C to +40°C

www.ashvision.com

At Ash we design, develop and manufacture all our user centric solutions in-house and are proud of our award winning innovation process. We use creative Design Thinking to actively empathise with our customers to understand their real unmet needs and jobs to be done. We seek meaningful engagement and co-creation with our end users so we can develop the best possible solutions and services in the quality assurance industry resulting in cost savings, increased workflow efficiency, waste reduction and an overall improved quality process.

Contact Us

Ireland

ASH
B5, M7 Business Park,
Naas, Co. Kildare,
W91 P684, Ireland.

P: + 353-45-882212

F: + 353-45-882214

United Kingdom

ASH
Covert Farm,
Long Lane,
East Haddon,
Northamptonshire
NN6 8DU
UK

P: +44 (0)7592 523 767

E: info@ashvision.com

W: www.ashvision.com

