

HAND-HELD MAGNIFIERS AND READING MAGNIFIERS

OPTIMUM VISION COMFORT
- ALWAYS READY TO HAND.



HAND-HELD MAGNIFIERS & READING MAGNIFIERS



reddot design award
winner 2011

ILLUMINATED POCKET MAGNIFIERS

mobilux® LED

- Batteries are easy to change as the battery cover is easy to open and contact springs clearly indicate the battery position.
- Improved comfort due to eyelet for lanyard.
- The most suitable colour temperature can be preselected via the 2 attachable filters supplied.
- The LED light does not require changing.
The LED light used by Eschenbach Optik has a service life of approx. 50,000 hours.
- The battery life is up to 10x longer in comparison to a conventional light.
- The voltage is controlled via a step-up converter, which produces consistent illumination even if the voltage drops due to low batteries.
- Bright, even illumination provided by SMD LED and collimation optics.
- Defined focusing distance: in the case of all aspheric magnifiers and reading magnifiers from Eschenbach Optik the focusing distance (dE: distance from the eye to the virtual image) is precisely 400 mm. The image can therefore be seen clearly with readers which have an addition of +2.5 dpt.

The mobilux® DIGITAL can be found on Page 26.




CUSTOMER BENEFITS

- ✚ Suitable for reading small print, e.g. in a telephone directory or on price tags.
- ✚ Simple to use – large, easy to operate light switch and eyelet for lanyard.
- ✚ Choice of 3 different lights colours via filters.
- ✚ The batteries are easy to change.
- ✚ With handy eyelet on the end of the handle for round-the-neck cord.
- ✚ Incl. protective case for lens.

Item	Dimensions	Magnification	dpt
15112	Ø 60 mm	3,0 x	12,0
15113	75 x 50 mm	3,5 x *	10,0
15114	75 x 50 mm	4,0 x **	16,0
151141	Ø 60 mm	4,0 x	16,0
15115	Ø 58 mm	5,0 x	20,0
15116	Ø 58 mm	6,0 x	24,0
15117	Ø 35 mm	7,0 x	28,0
151110	Ø 35 mm	10,0 x	38,0
151112	Ø 35 mm	12,5 x	50,0

* Commercial magnification

** Diffractive aspheric lens without cera-tec® coating

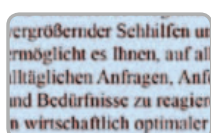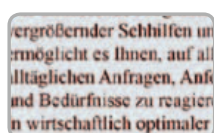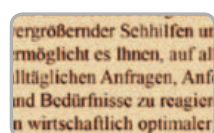

Accessories: Lanyard Order No. 15981



A GLOBAL BRAND WITH TRADITION

Since its establishment in 1913, Eschenbach Optik has been combining traditional business values with modern innovative flair. Within a century Eschenbach Optik has become one of the leading suppliers of optical products - both in Germany and abroad. The product range extends from eyewear frames and sunglasses, to binoculars and right through to magnifiers, telescopic glasses and special vision aids.

With subsidiaries in twelve countries, more than 30 international trade partners and over 630 employees throughout the world, today Eschenbach Optik has a global presence.

VISION TECHNOLOGY PRODUCTS

In this catalogue you will find our extensive range of magnifying vision aids. The term "Vision Technology Products" stands for innovative research and development, modern technologies, as well as precision production methods - and a level of quality which meets the very highest demands in terms of optimum optical performance, reliability and intuitive operation.

The diversity of our range of magnifying vision aids allows you to provide your customers with the optimum solution for meeting all their everyday vision requirements.

LOW VISION

In our catalogue you will therefore find a wide range of products to meet the highly diverse needs of those suffering from visual impairment, making their daily life so much easier and therefore enabling them to regain maximum independence and maximum quality of life.

EASY VISION

As we get older, our vision deteriorates and this is something that happens to everyone. From the age of 60, this change becomes increasingly noticeable and sometimes even the best corrective eyewear is not enough.

With the "Easy Vision" concept from Eschenbach Optik, additional, innovative optical products are able to compensate for the visual impairment and relaxed vision is once again possible.

MADE IN GERMANY

From the initial concept right through to the finished product - Eschenbach Optik carries out all of this work at its site in Nuremberg and this is how we are able to guarantee exceptional brand quality "MADE IN GERMANY".



SALES SUPPORT

On Pages 60-69, you can find our range of sales support products - including versatile displays and attractive presentation systems. Ask your representative about further advertising opportunities.

DETAILED INFORMATION

The "Technical Information" section contains a summary of technical specifications and you can find information on visual impairment in the "Frequently Asked Questions" section.

INDIVIDUAL LENS RANGE

Single vision or multifocal spectacle lenses with corrections to suit specific individual requirements, cut-off filters with individual strengths and reading glasses with individually adapted high additions can be found in the separate "2013 Individual Lens Range" catalogue.



Buy online at www.pepleroptics.com

PeplerOptics
Quality magnified.

Art.No
KAT-VTP03



 **ESCHENBACH**

Eschenbach Optik GmbH
Schopenhauerstrasse 10
90409 Nuremberg

Tel.	+49 (0) 911 3600-0
Fax	+49 (0) 911 3600-358
Internet	www.eschenbach-optik.com
E-Mail	mail@eschenbach-optik.com